


# TurboID and miniTurbo- engineered promiscuous biotin ligases for efficient proximity labeling in living cells and organisms

Researchers at Stanford have engineered two promiscuous biotin ligases for non-toxic, efficient proximity labeling (PL) in living cells and organisms. PL is a powerful technique for the proteomic analysis of macromolecular complexes, organelles or protein interaction networks. In PL, a promiscuous labeling enzyme is fused to a protein of interest or targeted to a subcellular region and upon the addition of biotin, the enzyme begins labeling nearby (within a few nanometers) endogenous proteins. Then the biotin labelled proteins are harvested and identified by mass spectrometry. Existing PL techniques are not optimal as they require either toxic reagents or long labelling times that prevent the study of dynamic structures or processes. To overcome these limitations the inventors have engineered two mutants of biotin ligase, called TurboID and miniTurbo, which are both more active and efficient than existing variants. TurboID and miniTurbo enable PL with greater signal in shorter time windows under physiologically relevant conditions in live cells. This technology will expand the scope and utility of PL by enabling the study of dynamic processes *in vivo* and thus can provide insight into proteomic changes to aid general research and therapeutic development.

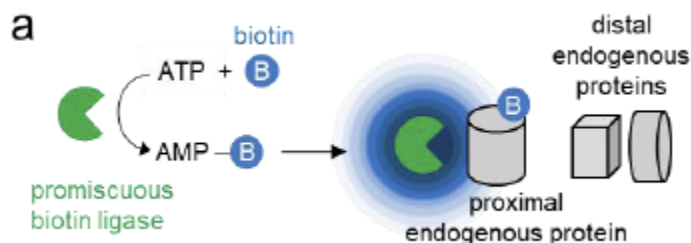

Proximity-dependent biotinylation by promiscuous biotin ligases. Ligases catalyze the formation of biotin-5'-AMP anhydride, which diffuses out of the active site to biotinylate proximal endogenous proteins.

### **Stage of research**

In addition to shortening PL time by 100-fold and increasing PL yield in cell culture, the inventors have shown that TurboID enabled biotin-based PL in new settings including yeast, *Drosophila* and *C. elegans*.

## **Applications**

- Proteomic analysis for:
  - Research
  - Therapeutic development

## **Advantages**

- Efficient- shortens PL time by 100-fold
  - Short labeling time-can be used to resolve dynamic proteomic changes during development or in response to drugs or stimuli
- Enhanced activity- proteomic labeling in just 10 mins (rather the 18-24hrs required by other methods)
- Non-toxic- can be used in living cells and organisms
  - Allows proteomic mapping under physiologically relevant conditions with high temporal resolution
- Simple labeling protocol
- Can be genetically targeted to any subcellular compartment or fused to any protein of interest
- Allows clinically relevant studies in animal models

## **Publications**

- Branon, T., Bosch, J., Sanchez, A. et al. [Efficient proximity labeling in living cells and organisms with TurboID](#). *Nat Biotechnol* 36, 880–887 (2018).

## Patents

- Published Application: [WO2019143529](#)
- Published Application: [20210214708](#)
- Issued: [12,110,522 \(USA\)](#)

## Innovators

- Alice Ting
- Tess Branon

## Licensing Contact

### **Mona Wan**

Senior Licensing Manager for Special Projects

[Email](#)